EN Operating Instructions

Thank you for purchasing WI-911 oral irrigator from B.Well. The appliance is helpful after meals, in case of tooth decay and plaque removal. It precludes dental tartar occurrence and is perfect for prevention of caries and periodontal diseases. WI-911 oral irrigator is a small unit for oral care if you have dental bridges, implants, crowns and brackets.

It can also be used for gum massage that improves blood circulation and promotes gum healing. Three operating modes ensure the maximum efficiency of the unit.

1. SAFETY PRECAUTIONS

① WARNING!

Never use the following as solutions for oral irrigation:

- Substances and solutions containing suspended particles (decoctions, suspensions, infusions, etc.). Such particles may cause obstruction of the suction hose and the nozzle. Oral irrigator solution is an organic solvent and its use may have a strong effect on the suction hose and components of the appliance.
- Antiseptic solutions (potassium permanganate, miramistin and other antiseptic drugs).
- Salt solutions, toothpaste or chemical substances.

Allowed for oral irrigation:

- Clean lake water (up to 40 °C). Use filtered water if you have poor quality tap water in your region.
- Special rinses for use with irrigator / common mouth rinse solution (always mixed in a definite ratio with water. For recommended proportions refer to the rinsing manufacturer's instructions).
- ① CAUTION! If special solutions for irrigator / common mouth rinse solution are applied, be sure to clean the device with fresh lake water after using such rinses in the water tank with water and run an extra irrigation cycle, by draining the appliance through the nozzle. Pay close attention to the expiry date of permitted solutions. Never use outdated solutions.

① CAUTION!

- The appliance has built-in batteries. Avoid throwing it into fire, heating, charging, using or leaving the unit in locations with high air temperature.
- Never use the adapter with damaged cord or power plug.

④ PRECAUTIONS:

- Do not use the appliance with the charging adapter in the bath or in the shower. Never leave, nor keep the charging adapter where it can fall or get into a bathtub, a sink or a WC pan. This may result in electrocution or fire.
- Do not try to reach an irrigator that has fallen into water; unplug first.
- Do not submerge the irrigator in water.
- Never use the appliance with damaged charging adapter, or if the power cord is loosely connected to the power socket. This may result in electrocution, short-circuit or fire.
- Never break, bend, pull or twist the power cord with excessive force. Do not place heavy objects onto the power cord and do not place the cord between such objects. This may result in electrocution or fire. Do not submerge the charging adapter into water and never wash it with water. This may result in electrocution or short-circuit.
- Never plug in / disconnect the adapter from the power socket with wet hands.
- Store the appliance out of reach of children. The failure to observe the rule may result in accidents and hazardous situations, for instance, accidental swallowing of the accessories or detachable parts of the unit by a child.
- During maintenance, plug the charging device out of the power socket. The failure to observe this may result in electrocution.
- Regularly clean the power plug and the charging adapter to avoid dust ingress on them. The failure to observe the rule may result in fire due to damaged insulation caused by humidity and dust accumulation.
- At all times, make sure that the device works from a power source with the relevant rated voltage. The failure to observe the rule may result in electrocution or fire.
- Never modify the appliance, disassemble or repair it on your own. This may result in fire, electrocution or an injury. If repair is needed (for instance, to replace a battery) contact the authorised service centre of B.Well.
- Do not store the unit where it can fall and break into pieces (e.g. on a sink).
- Do not use the appliance with lost nozzles.
- Do not use the appliance connected to the power socket.

④ WARNING:

- Children and persons who are unable to operate the device, and people with weak oral cavity sensibility should not use the appliance. The failure to observe the rule may result in electrocution or fire.
- Never leave the unit in the water tank with temperature above 40°C.
- Always use the correct charging adapter from the power socket, holding the former onto the plug, rather than the cord. The failure to observe the rule may result in electrocution, short-circuit or fire.
- Keep the appliance clean, use only fresh and lake water for mouth cleaning. Some types of liquids may damage to the appliance. Therefore, never use liquids from the forbidden list for irrigation.
- After use never forget to drain the remaining water from the tank and operate the appliance for 2-3 seconds until no more water appears. If no health and hygiene rules are observed, this may result in growth of fungi and bacteria, which are harmful for health. This may also cause odour and malfunctions to originate.
- The appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, or persons without relevant expertise, unless they are supervised by a person responsible for their safety, and provided they have been given precise use instructions. Be sure the appliance is not within reach of children for play.

2. INTENDED USE

- WI-911 oral irrigator with accessories is intended to prevent:
 - gum bleeding, caries, periodontal diseases, gingivitis, periimplant pathology; medicated irrigation and massaging of the oral mucosa;
 - for effective cleaning of orthodontic systems, and for care of oral cavity and space between teeth.
 - it is intended for private use.

3. SCOPE OF SUPPLY

- Oral irrigator
- 2 standard nozzles
- Water tank
- Charger (adapter)
- 5 Charging socket caps (2 pcs)
- 6 Battery kit (a cell)
- Operating instructions
- Warranty card
- Box

- ① Note: There may be moisture inside the water tank or suction hose at the time of purchase. This is moisture left from the distilled water used for product testing and is safe for humans.

STANDARD NOZZLES FOR WI-911 irrigator

Clean space between teeth from bacteria and food particles. Effectively massages gums. Suitable for cleaning bridges, implants, crowns and brackets.

Recommended replacement interval is 6 months.

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5. CHARGING THE APPLIANCE

- Plug the adapter cord into the irrigator charging socket and connect the device to a power socket (make sure no moisture is inside the irrigator socket).
- CH indicator will flash while the irrigator is charging; once the irrigator is charged, CH indicator will turn off.
- Close the charging socket with a cap after the charging is complete.

- ① Note: If CH indicator keeps flashing during and after the procedure, it means the irrigator battery needs to be charged. Battery charging is completed in some 8 hours.

6. USING THE APPLIANCE


- Attach the nozzle into the relevant recess of the main unit. To extract the nozzle, press the nozzle release button and pull the nozzle upwards.
- Fill in the tank with water by any one of the suggested methods:

- Open the water tank cap, hold the tank with water at temperature of 40°C or less. Close the cap making sure it is tightly closed.
- Reopen the water tank, while holding the main unit upright. Pour the tank with water at temperature of 40°C or less and reinsert the tank into the appliance.

- ① CAUTION! Never use a water-filled irrigator without the nozzle, which may cause splashing.

- Use only fresh water as all times.
- Never pour solutions from the forbidden list, specified in the beginning of the instruction. See 1. Safety Precautions.
- To ensure trouble-free operation make sure the irrigator is turned off when the tank is being filled with water.
- Do not operate the appliance after use. Operating the appliance without water may result in a malfunction.

- Before pressing the power switch place the nozzle into your mouth. Bend a water tank, and hold the appliance upright at the same direction the nozzle onto the teeth. Slightly close your mouth during irrigation to prevent splashing.

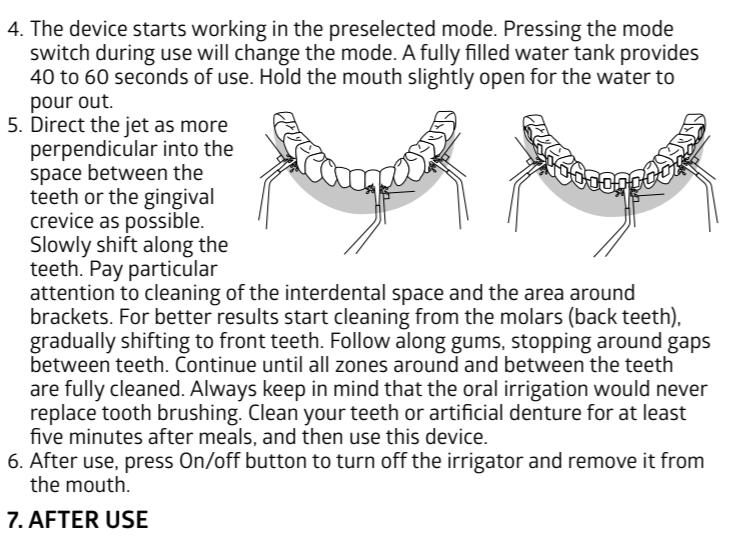4. The device starts working in the preselected mode. Pressing the mode switch during use will change the mode. A fully filled water tank provides 40 to 60 seconds of use. Hold the mouth slightly open for the water to pour out.

5. Direct the jet as more perpendicular to the teeth or the gingival crevice as possible.

Slowly shift along the teeth. Pay particular attention to cleaning of the interdental space and the area around brackets. For better results start cleaning from the molars (back teeth), gradually shifting to front teeth. Follow along gums, stopping around teeth with fillings. Clean your teeth on artificial denture for at least five minutes after meals, and then use this device.

6. After use, press On/off button to turn off the irrigator and remove it from the mouth.

7. AFTER USE

- Without detaching the irrigator nozzle, open the water tank cap or detach the tank and drain water.
- Turn on the appliance for a few seconds to drain the remaining water.
- Wipe the irrigator with a dry cloth.

① NOTE:

- Always drain the remaining water. This will help prevent contamination and growth of bacteria in water that remained in the device.

- If you intend not to use the device for a long time, do not forget to wipe out and dry the appliance before storage.

8. CARE AND CLEANING

Cleaning Turn off the device, detach the water tank and the nozzle. When washing, never use water with temperature above 50°C.

Use only water or PH neutral detergents to clean the device. Never use caustic soda or abrasives for cleaning since those can damage it.

Main Unit

Use a cloth to clean the main unit. Use liquid soap or water to remove hard-to-clean dirt. Use only PH neutral detergents. Never use caustic and abrasive substances (e.g. vinegar or lime scale solvents), since it may do harm to the appliance.

Never immerse the main unit into water.

After washing, wipe the main unit dry with a cloth.

Water tank

Rinse the tank with water.

If you do not intend to use the appliance for a week or longer, dry the tank from the inside and outside.

Do not submerge the tank in water.

Use only soft PH neutral detergent for cleaning.

Nozzle

Wipe with water and wipe out with a soft cloth.

The recommended nozzle replacement interval is 6 months.

9. SUCȚIAZĂ

• Rinse with water, wipe out with a soft cloth.

• Do not bend, pull or twist.

• Rinse with water and wipe out with a soft cloth.

• The recommended nozzle replacement interval is 6 months.

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• The irrigator does not work. You have not inserted the battery or it has not been used for more than 3 months.

• The irrigator is not charged. The battery duty cycle is about to end.

• The battery is not sufficiently charged.

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